

White Label PCR Consumables

The White Label PCR Consumables are available as regular profile tubes of 0.2mL volume capacity, which is a great choice for standard PCR/qPCR thermocyclers, strips of 8-tubes with caps and 96-well plates.

All White Label PCR Consumables are made out of polypropylene, with a specially polished interior which significantly prevents any protein adhesion on the tube walls, enabling you to fully recover your sample.

Certified DNA, DNase, RNA, RNase and Pyrogen-free, White Label PCR Consumables are the perfect choice for the wide majority of applications.

Thanks to the high-quality polypropylene (PP) material from which they are constructed, the White Label PCR Consumables can also be autoclaved, if required.

White Label PCR Consumables for the Ultimate Heat Distribution

Due to their thin-walled construction, the White Label PCR Consumables offer the ultimate transmission and distribution of heat. Designed to be universally compatible with the majority of the thermocyclers on the market, White Label PCR Consumables are the optimal choice for worldwide laboratory set-ups.

The V-shape bottom of the White Label PCR Consumables enables full sample recovery and, along with the thin walls and leak-proof closure, the tube design minimises any condensation. The evaporation rate of the White Label PCR Tubes is less than 3%.


In addition to these features, the ergonomic flat cap design enables one-hand opening and closing, increasing your work efficiency, comfort and throughput.

Highlights:

- Thin-walled
- Even temperature transmission
- Minimal condensation
- Evaporation less than 3%
- Universal thermocycler compatibility
- Low and regular profile tubes
- One-hand operation
- Certified DNA, DNase, RNA, RNase and Pyrogen-free
- Autoclavable

Ordering Information

Cat. No.	Description	Dimensions	Packaging
WL3010800	PCR tube 0,2mL, clear, flat cap	0,2mL	10x1000 pcs.
WL3037020	PCR strips of 8 tubes 0,2mL, attached flat cap	8x0,2mL	5x10x12 pcs.
WL3047010	(q)PCR plate 96-well, 0,2mL, non-skirted, clear	96x0,2mL	5x5x10 pcs.
WL3048020	(q)PCR plate sealing membrane (adhesive film)		10x10 pcs.

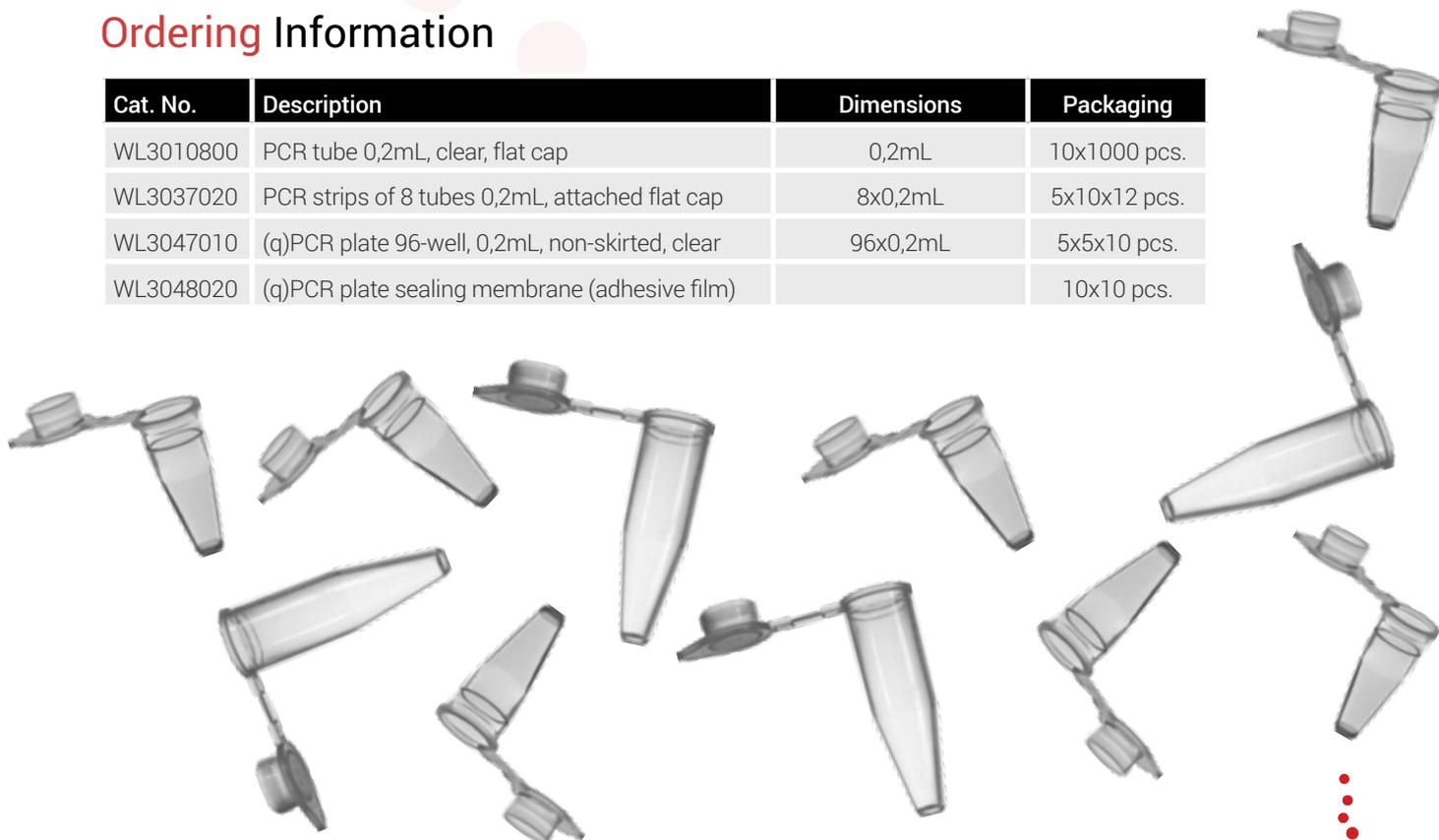

White Label Microtubes

The White Label Microtubes are made out of a transparent polypropylene of the highest clarity for the best visibility of your sample.

Manufactured by a top-class production facility and fully compliant with international norms and standards, these microcentrifuge tubes have highly-polished interiors, keeping the risk of any protein adhesion to an absolute minimum.

Designed to enable one-hand operation when opening and closing, as well as the ultimate sample sealing and evaporation protection, these White Label Microtubes are the optimal choice for the laboratories worldwide.

Ultimate Work Comfort with White Label Microtubes

The conical-base-design and easy-to-read graded measurement markings range from 0.1 to 1.5mL, providing you with the ultimate work comfort when freezing your samples.

The microtubes are equipped with a frosted hinged lid which makes the labelling of your samples extremely easy.

Highlights:

- Frosted hinged lid for easy labelling
- Easy-to-read measurement markings
- Single-handed operation
- Conical base
- Perfect sealing and evaporation protection
- Polished interior prevents protein adhesion

Ordering Information

Cat. No.	Description	Dimensions	Packaging
WL3116C50	Microcentrifuge tube, 0,5mL, sterile, clear	0,5mL	5x1000 pcs.
WL3205800	Microcentrifuge tube, 1,5mL, clear	1,5mL	16x500 pcs.
WL3204800	Microcentrifuge tube, 2,0mL, clear	2,0mL	16x500 pcs.



White Label Pipette Tips

Manufactured at a top-class production facility, compliant with international norms and standards, and having more than 20 years of experience in liquid handling consumables, the White Label Pipette Tips offer the ultimate product portfolio for worldwide laboratories.

Universally Compatible, Certified, Low-Retention Pipette Tips

Made of a high-quality virgin polypropylene and produced on top-class Swiss-made injection moulds, White Label Pipette Tips offer a universal compatibility with the majority of pipette brands on the market.

Designed to sustain a great dimensional consistency, having a great finish, flash-free orifice and crystal clarity, the White Label Pipette Tips only retain the absolute minimum of your sample.

All White Label Pipette Tips are certified RNase, DNase and Pyrogen-free.

Available in a large range of racked and bulk packaging options as well as filter and non-filter tips, the White Label Pipette Tips Solution will cover the wide scope of demand worldwide.


Ideal Choice for Molecular Biology, Clinical and Life Science

Equipped with filters made of pure HDPE (High Density Polyethylene) and free from any cellulose additives, these White Label Filter Tips effectively prevent cross-contamination on critical applications.

Having the filter pore size three times smaller than the common market standard, the White Label Filter Tips are the optimal choice for Molecular Biology as well as Clinical and Life Science.

Highlights:

- RNase, DNase and Pyrogen-free
- Made of virgin polypropylene
- Available as filter tips
- Filter pore size is three times smaller than the market average
- Produced on Swiss-made injection moulds.
- Universally compatible with most pipette brands on the market
- Great finish
- Flash-free orifice
- Crystal clarity
- Great dimensional consistency

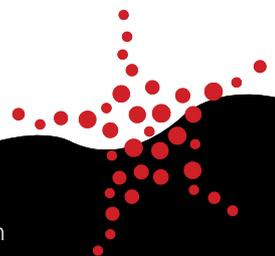


Ordering Information

Pipette tips, non-filtered, bulk and racked

Cat. No.	Description	Dimensions	Packaging
WLB1002810	Pipette tip 0.1-10 µL, clear	0.1-10 µL	20x1000 pcs.
WLB1002C50	Pipette tip 0.1-10 µL, clear, racked	0.1-10 µL	5x10x96 pcs.
WLB1003C50	Pipette tip 0.1-10 µL, clear, racked, sterile	0.1-10 µL	5x10x96 pcs.
WLB1062810	Pipette tip 0.1-10µL XL/20µL clear	0.1-10 µL	20x1000 pcs.
WLB1062C50	Pipette tip 0.1-10 µL XL/20µL clear, racked	0.1-10 µL	5x10x96 pcs.
WLB1063C50	Pipette tip 0.1-10 µL XL/20µL, clear, racked, sterile	0.1-10 µL	5x10x96 pcs.
WLB1121860	Pipette tip 1-200 µL, clear	1-200 µL	30x1000 pcs.
WLB1121861	Pipette tip 1-200 µL, yellow	1-200 µL	30x1000 pcs.
WLB1121C50	Pipette tip 1-200 µL, clear, racked	1-200 µL	5x10x96 pcs.
WLB1121C51	Pipette tip 1-200 µL, yellow, racked	1-200 µL	5x10x96 pcs.
WLB1122C50	Pipette tip 1-200 µL, clear, racked, sterile	1-200 µL	5x10x96 pcs.
WLB1122C51	Pipette tip 1-200 µL, yellow, racked, sterile	1-200 µL	5x10x96 pcs.
WLB1127800	Pipette tip 1-300 µL, clear	1-300 µL	10x1000 pcs.
WLB1127C50	Pipette tip 1-300 µL, clear, racked	1-300 µL	5x10x96 pcs.
WLB1128C50	Pipette tip 1-300 µL, clear, racked, sterile	1-300 µL	5x10x96 pcs.
WLB1202800	Pipette tip 100-1000 µL, clear	100-1000 µL	10x1000 pcs.
WLB1202802	Pipette tip 100-1000 µL, blue	100-1000 µL	10x1000 pcs.
WLB1202C40	Pipette tip 100-1000 µL, clear, racked	100-1000 µL	4x8x96 pcs.
WLB1202C42	Pipette tip 100-1000 µL, blue, racked	100-1000 µL	4x8x96 pcs.
WLB1203C40	Pipette tip 100-1000 µL, clear, racked, sterile	100-1000 µL	4x8x96 pcs.
WLB1203C42	Pipette tip 100-1000 µL, blue, racked, sterile	100-1000 µL	4x8x96 pcs.
WLB1204800	Pipette tip 100-1250 µL XL, clear	100-1250 µL XL	10x1000 pcs.
WLB1204C40	Pipette tip 100-1250 µL XL, clear, racked	100-1250 µL XL	4x8x96 pcs.
WLB1205C40	Pipette tip 100-1250 µL XL, clear, racked, sterile	100-1250 µL XL	4x8x96 pcs.

Pipette tips 5mL and 10mL

Cat. No.	Description	Size	Packaging
WLB1218610	Pipette tip 1-5 mL, clear (Capp, Gilson, HTL)	1-5 mL	10x250 pcs.
WLB1222610	Pipette tip 1-5 mL, clear, sterile (Capp, Gilson, HTL)	1-5 mL	10x250 pcs.
WLB1220500	Pipette tip 1-10 mL, clear (Capp, Gilson, Eppendorf)	1-10 mL	5x200 pcs.
WLB1224500	Pipette tip 1-10 mL, clear, sterile (Capp, Gilson, Eppendorf)	1-10 mL	5x200 pcs.

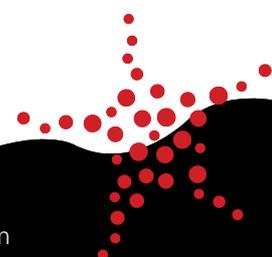

Filter tips, bulk and racked


Cat. No.	Description	Size	Packaging
WLB2002810	Pipette Filter Tips 10 µL, clear, bag	10 µL	20x1.000 pcs.
WLB2003C50	Pipette Filter Tips 10 µL, clear, sterile, racked	10 µL	5x10x96 pcs.
WLB2062810	Pipette Filter Tips 10 µL XL/ 20 µL, clear, bag	10 µL XL/ 20 µL	20x1.000 pcs.
WLB2063C50	Pipette Filter Tips 10 µL XL/ 20 µL, clear, sterile, racked	10 µL XL/ 20 µL	5x10x96 pcs.
WLB2117800	Pipette Filter Tips 20 µL, clear, bag	20 µL	10x1000 pcs.
WLB2118C50	Pipette Filter Tips 20 µL, clear, sterile, racked	20 µL	5x10x96 pcs.
WLB2119800	Pipette Filter Tips 100 µL, clear, bag	100 µL	10x1000 pcs.
WLB2120C50	Pipette Filter Tips 100 µL, clear, sterile, racked	100 µL	5x10x96 pcs.
WLB2127800	Pipette Filter Tips 200 µL, clear, bag	200 µL	10x1000 pcs.
WLB2128C50	Pipette Filter Tips 200 µL, clear, sterile, racked	200 µL	5x10x96 pcs.
WLB2134800	Pipette Filter Tips 300 µL, clear, bag	300 µL	10x1000 pcs.
WLB2135C50	Pipette Filter Tips 300 µL, clear, sterile, racked	300 µL	5x10x96 pcs.
WLB2202800	Pipette Filter Tips 1000 µL, clear, bag	1000 µL	10x1000 pcs.
WLB2203C40	Pipette Filter Tips 1000 µL, clear, sterile	1000 µL	4x8x96 pcs.
WLB2204800	Pipette Filter Tips 1250 µL, clear, bag	1250 µL	10x1000 pcs.
WLB2205C40	Pipette Filter Tips 1250 µL, clear, sterile, racked	1250 µL	4x8x96 pcs.

Reload tips

Cat. No.	Description	Size	Packaging
WLB1002C1R	Pipette Tips 10 µL, clear, reload	10 µL	10x10x96 pcs.
WLB1062C1R	Pipette Tips 10 µL XL/ 20 µL, clear, reload	10 µL XL/ 20 µL	10x10x96 pcs.
WLB1121C1R	Pipette Tips 200 µL, clear, reload	200 µL	10x10x96 pcs.
WLB1127C2R	Pipette Tips 300 µL, clear, reload	300 µL	10x8x96 pcs.
WLB1202C3R	Pipette Tips 1000 µL, clear, reload	1000 µL	10x5x96 pcs.

Empty tip racks

Cat. No.	Description	Size	Packaging
WLB4500C5E	Empty Rack Boxes 10µL /10 µL XL, clear black	10 µL/10 µL XL	5x10 pcs.
WLB4501C5E	Empty Rack Boxes 200 µL /300 µL, clear black/yellow	200 µL/300 µL	5x10 pcs.
WLB4503C4E	Empty Rack Boxes 1000 µL, clear black/blue	1000 µL	4x8 pcs.
WLB4504C4E	Empty Rack Boxes 1250 µL, clear black/blue	1250 µL	4x8 pcs.

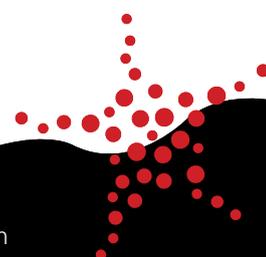


Pipette tips, non-filtered, bulk and racked, low retention

Cat. No.	Description	Size	Packaging
WLB4002810	Low Retention Tips 10 µL, clear, bag	10 µL	20x1.000 pcs.
WLB4003C50	Low Retention Tips 10 µL, clear, sterile, racked	10 µL	5x10x96 pcs.
WLB4062810	Low Retention Tips 10 µL XL, clear, bag	10 µL XL	20x1.000 pcs.
WLB4063C50	Low Retention Tips 10 µL XL, clear, sterile, rack	10 µL XL	5x10x96 pcs.
WLB4121860	Low Retention Tips 200 µL, clear, bag	200 µL	10x1000 pcs.
WLB4122C50	Low Retention Tips 200 µL, clear, sterile, racked	200 µL	5x10x96 pcs.
WLB4127800	Low Retention Tips 300 µL, clear, bag	300 µL	10x1000 pcs.
WLB4128C50	Low Retention Tips 300 µL, clear, sterile, racked	300 µL	5x10x96 pcs.
WLB4202800	Low Retention Tips 1000 µL, clear, bag	1000 µL	10x1000 pcs.
WLB4203C40	Low Retention Tips 1000 µL, clear, sterile, racked	1000 µL	4x8x96 pcs.
WLB4204800	Low Retention Tips 1250 µL, clear, bag	1250 µL	10x1000 pcs.
WLB4205C40	Low Retention Tips 1250 µL, clear, sterile, racked	1250 µL	4x8x96 pcs.

Filter tips, bulk and racked, low retention

Cat. No.	Description	Size	Packaging
WLB5003C50	Low Retention Filter Tips 10 µL, clear, sterile	10 µL	5x10x96 pcs.
WLB5063C50	Low Retention Filter Tips 10 µL XL, clear, sterile	10 µL XL	5x10x96 pcs.
WLB5118C50	Low Retention Filter Tips 20 µL, clear, sterile	20 µL	5x10x96 pcs.
WLB5120C50	Low Retention Filter Tips 100 µL, clear, sterile	100 µL	5x10x96 pcs.
WLB5128C50	Low Retention Filter Tips 200 µL, clear, sterile	200 µL	5x10x96 pcs.
WLB5135C50	Low Retention Filter Tips 300 µL, clear, sterile	300 µL	5x10x96 pcs.
WLB5203C40	Low Retention Filter Tips 1000 µL, clear, sterile	1000 µL	4x8x96 pcs.
WLB5205C40	Low Retention Filter Tips 1250 µL, clear, sterile	1250 µL	4x8x96 pcs.

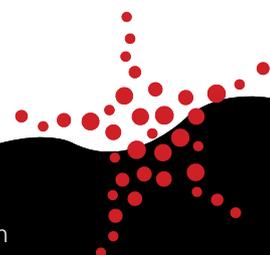

White Label Reagent Reservoirs

The White Label Reagent Reservoirs are made of a unique blend of Polystyrene and provide a sturdy, universal reservoir solution to any global laboratory.

Large Sturdy Base with Graduation Marks

Available in three different sizes - 10mL, 25mL, 50mL - the White Label Reagent Reservoirs provide a large, sturdy base with the broad edges, standing firmly on your workbench.

The graduation marks enable accurate dispensing and sample volume identification

Angled Interior and Spillproof

The angled interior of the White Label Reagent Reservoirs supports the collection of liquid and decreases reagent waste.

While processing liquid there is no need to worry about spillage because the reservoirs are spillproof.

Compact Storage and Easy Pouring

To store any remaining liquid, just empty the reservoirs using one of the corners.

For compact storage you can stack the reservoirs easily and separate them by using the reservoir lip.

These White Label Reagent Reservoirs are ideal for using with all Single and Multichannel Pipettes.

To suit the various pipettes, the White Label Reagent Reservoirs come in three different sizes:

- 10mL, V-shaped bottom (ideal with 8- and 16-channel pipettes)
- 25mL, flat bottom (ideal with multichannel pipettes)
- 50mL, V-shaped bottom (ideal with 12-channel pipettes)

Highlights:

- Multiple volume options
- Large, sturdy base
- Spillproof
- Perfect for multichannel pipettes
- Compact storage



Ordering Information

Cat. No.	Description	Dimensions	Packaging
WL3701119	Reagent Reservoirs 10mL, bulk	10mL	5x100 pcs.
WL3702119	Reagent Reservoirs 25mL, bulk	25mL	5x100 pcs.
WL3703119	Reagent Reservoirs 50mL, bulk	50mL	5x100 pcs.